



The total face experience

Introducing the FitLife mask

FitLife is a total face mask that doesn't contact the areas on and around the nose. FitLife seals around the perimeter of the face, making it a good solution for nasal area sensitive patients, hard-to-fit patients and claustrophobic patients.

For the medical professional

- With three sizes designed to help fit even the most challenging patient faces, FitLife helps you to find the right fit for your patient.
- Its over-molded cushion integrated to the mask faceplate and its three sizes makes your inventory management easy.
- FitLife helps to establish an effective seal around the patient's face, even under high pressure, making it suitable for a broad range of therapy.

For the patient

- FitLife was designed to reduce facial marks, potential skin sores and stress on the nose bridge, helping to improve therapy adherence.
- The full vision screen and lightweight body allows maximum visibility and increased comfort.
- The inlaid headgear tabs and snap clips make adjustments easier, for an easy set-up and a quick removal.

PHILIPS
RESPIRONICS

The difference is clear to see

FitLife premium low-stretch headgear for improved stability and a consistent fit

Large surface area equalizes pressure and flow inside the mask, minimizing eye irritation

Large inlaid headgear tabs for quick, easy adjustment and removal of headgear

Snap clips simplify set-up and eliminate the need to refit after removal

Entrainment valve allows access to room air if pressure drops below 3 cmH₂O



Increased performance and comfort


- Effortless and effective seal, even under high pressures
- Wide, full-vision screen for maximum visibility
- Optimum pressure distribution, minimizing pressure sores and skin breakdown

Simplicity

- Limited mask parts, due to over-moulded cushion integrated into the faceplate
- Three sizes including small for pediatric use (>7 years, 20 kgs)
- Easy-to-fit premium headgear with integrated inlaid tabs and snap clips

Putting FitLife to the test

In 2011, five sleep labs were asked to trial the FitLife mask*. The trial demonstrated that the majority of patients are likely to continue using FitLife for the duration of their treatment.

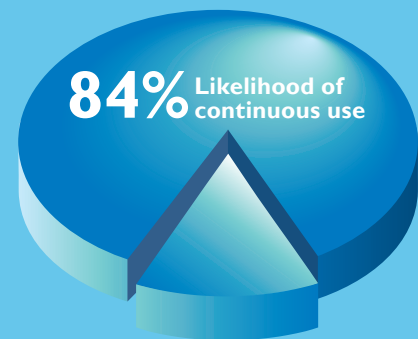

* Philips patient evaluation trial, Germany
FitLife vs leading competitor full face masks

FitLife description	Small	Large	X- large
Mask with headgear, EE version	1060803	1060804	1089994
Mask with headgear, SE version	1061711	1061712	1089995

Small size intended for both adult and pediatric use (7 years old or >20 kgs)

Please visit www.philips.com/respironics



© 2013 Koninklijke Philips N.V.
All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Broudy AE 7/4/13 MCI 4105690 PN 1063921

Philips Healthcare is part of Royal Philips

Philips Respironics Asia Pacific
+65 6882 5282
Philips Respironics Australia
+61 (2) 9947 0440
Philips Respironics United Kingdom
+44 800 1300 845

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668 USA
Customer Service
+1 724 387 4000
www.philips.com/respironics